


# 金银花-木犀草苷 (中国药典)

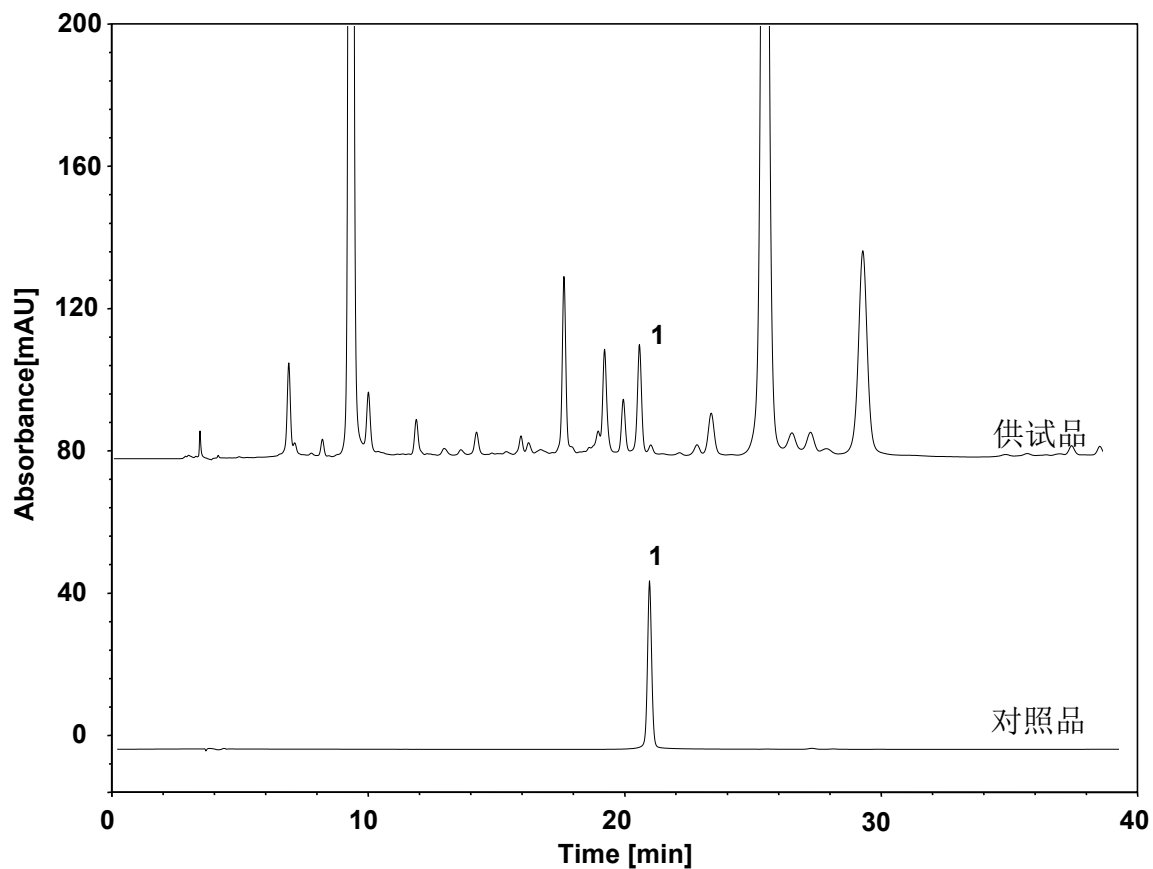

Column: ChromCore Phenyl, 5  $\mu$ m  
Dimension: 4.6  $\times$  250 mm  
Mobile Phase: A) 0.5%冰醋酸溶液  
B) 乙腈

Gradient:	t (min)	%A	%B
	0	90	10
	15	80	20
	30	80	20
	40	70	30

Flow Rate: 1.0 mL/min  
Temperature: 25  $^{\circ}$ C  
Injection: 10  $\mu$ L  
Detection: UV 350 nm  
Sample: 金银花  
Peak: 1. 木犀草苷