


# 羌活-特征图谱 (中国药典)

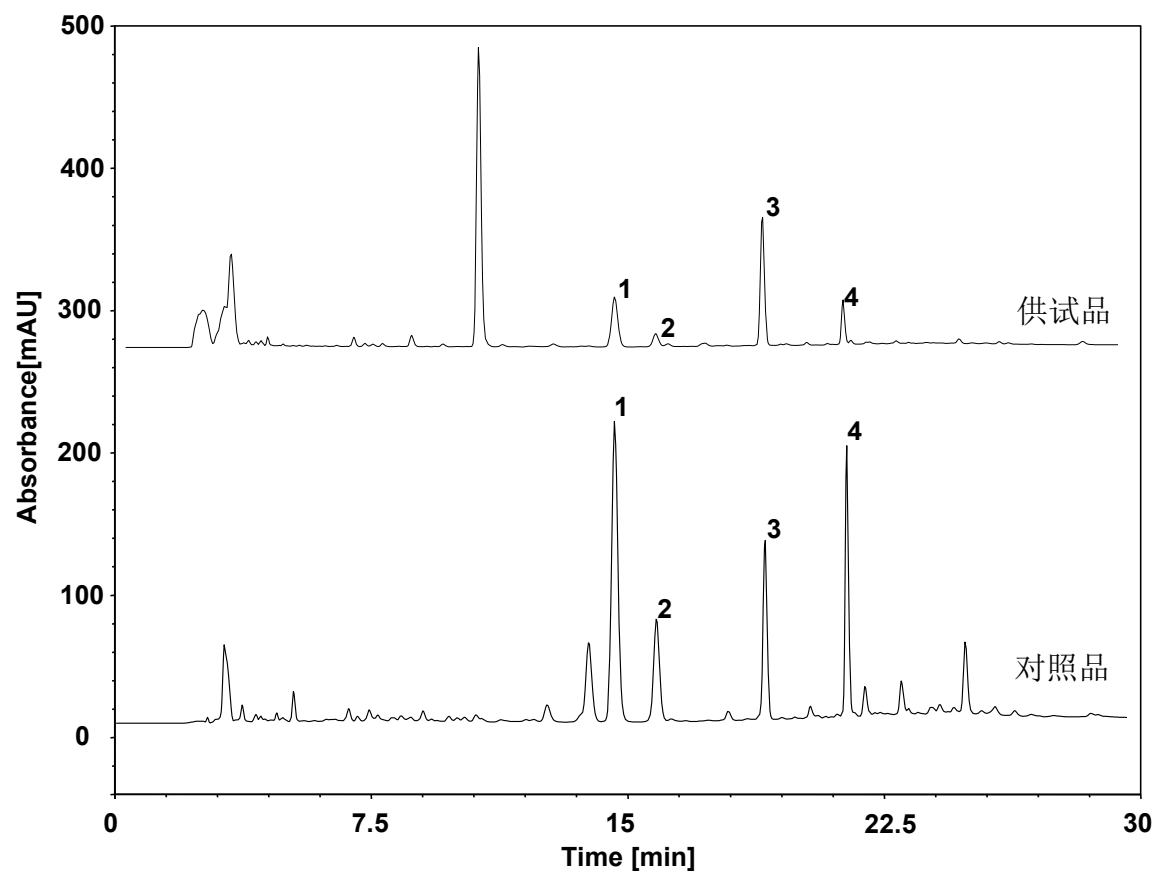

Column: ChromCore C18, 5  $\mu$ m

Dimension: 4.6  $\times$  250 mm

Mobile Phase: A) 乙腈

B) 0.1%磷酸溶液

Gradient:	t (min)	%A	%B
	0	48	52
	6	53	47
	12	53	47
	20	80	20
	30	80	20

Flow Rate: 1.0 mL/min

Temperature: 25  $^{\circ}$ C

Injection: 10  $\mu$ L

Detection: UV 246 nm

Sample: 羌活

Peaks:

1. 羌活醇

2. 阿魏酸苯乙醇酯

3. 异欧前胡素

4. 镰叶芹二醇