


奥美沙坦酯-溶出度

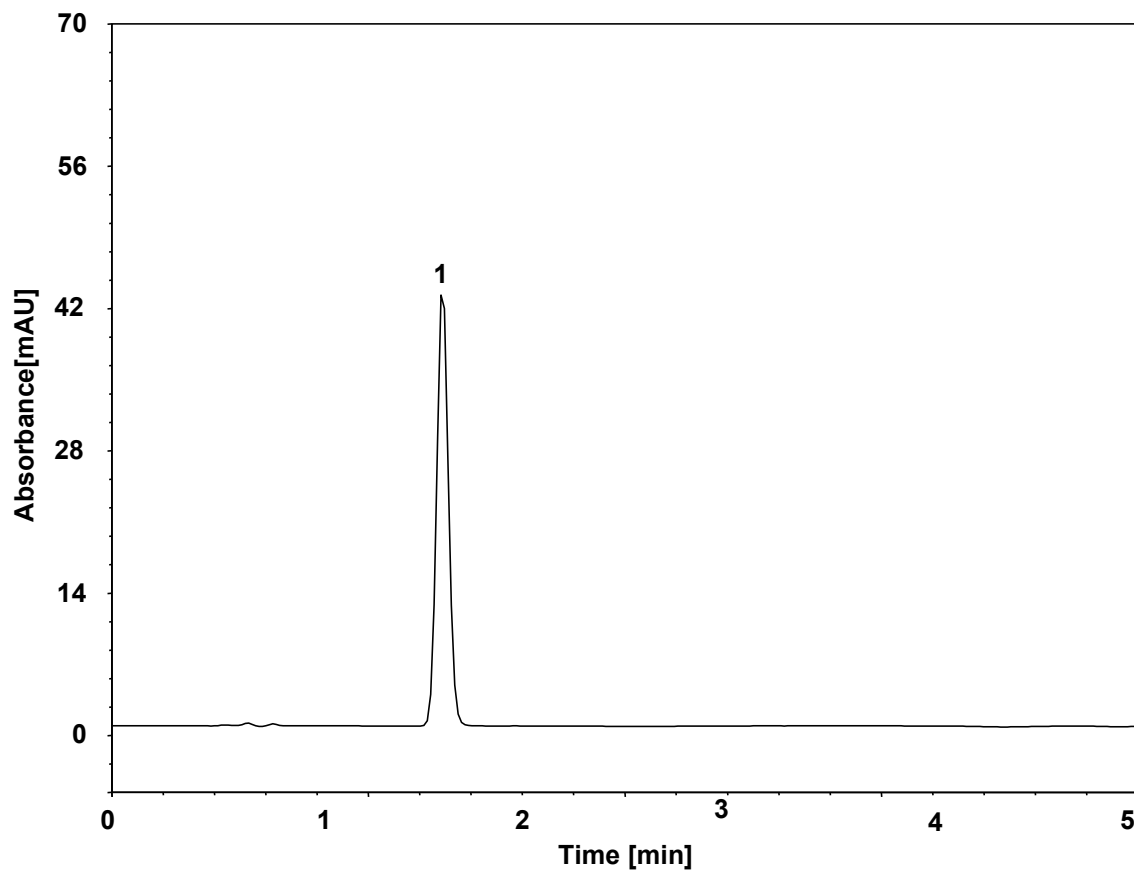

Column: ChromCore C18, 5 μ m
Dimension: 4.6 \times 100 mm
Mobile Phase: 40/60 v/v 乙腈/30 mM磷酸二氢钾溶液, pH2.5
Flow Rate: 2 mL/min
Temperature: 30 $^{\circ}$ C
Injection: 10 μ L
Detection: UV 254 nm
Sample: 奥美沙坦酯溶液
Peak: 1. 奥美沙坦酯

	保留时间 (min)	拖尾因子	理论塔板数
奥美沙坦酯	1.623	1.069	3545