


# 苯磺酸左氨氯地平

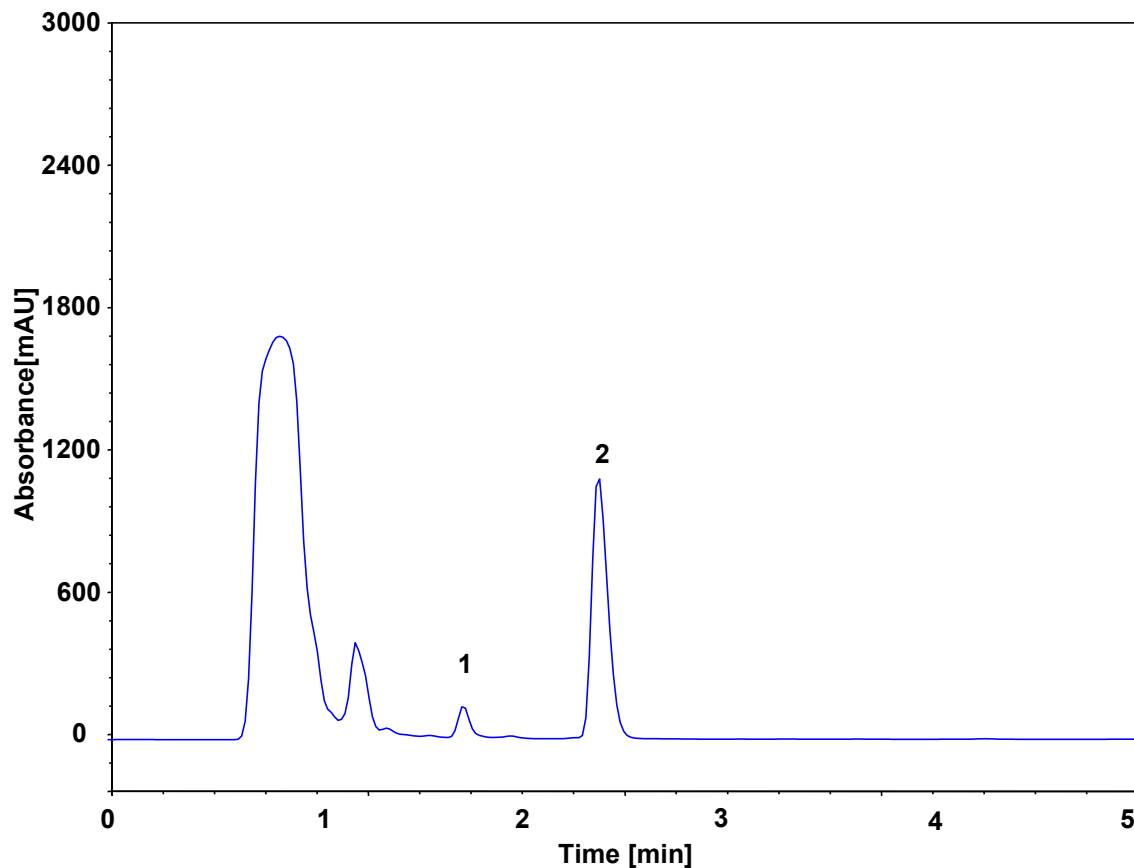

Column: ChromCore C18, 5  $\mu$ m  
Dimension: 4.6  $\times$  100 mm  
Mobile Phase: 40/60 v/v 乙腈/30 mM磷酸二氢钾溶液, pH2.5  
Flow Rate: 1.5 mL/min  
Temperature: 30  $^{\circ}$ C  
Injection: 50  $\mu$ L  
Detection: UV 237 nm  
Sample: 苯磺酸左氨氯地平系统适用性溶液  
Peaks: 1. 杂质  
2. 氨氯地平

|          | 保留时间  | 拖尾因子  | 分离度  | 理论塔板数 |
|----------|-------|-------|------|-------|
| 苯磺酸左氨氯地平 | 2.394 | 1.296 | 5.13 | 4284  |