


维生素B6



Column: ChromCore 120 C18, 5 μ m
Dimension: 4.6 \times 250 mm
Mobile Phase: 10/90 v/v 甲醇/0.04% 戊烷磺酸钠水溶液, pH3.0
Flow Rate: 1.0 mL/min
Temperature: 30 $^{\circ}$ C
Injection: 10 μ L
Detection: UV 291 nm
Peaks:
1. Vitamin B6 (维生素B6)
2. Vitamin B6 Impurity A (维生素B6杂质A)
3. Vitamin B6 Impurity B (维生素B6杂质B)

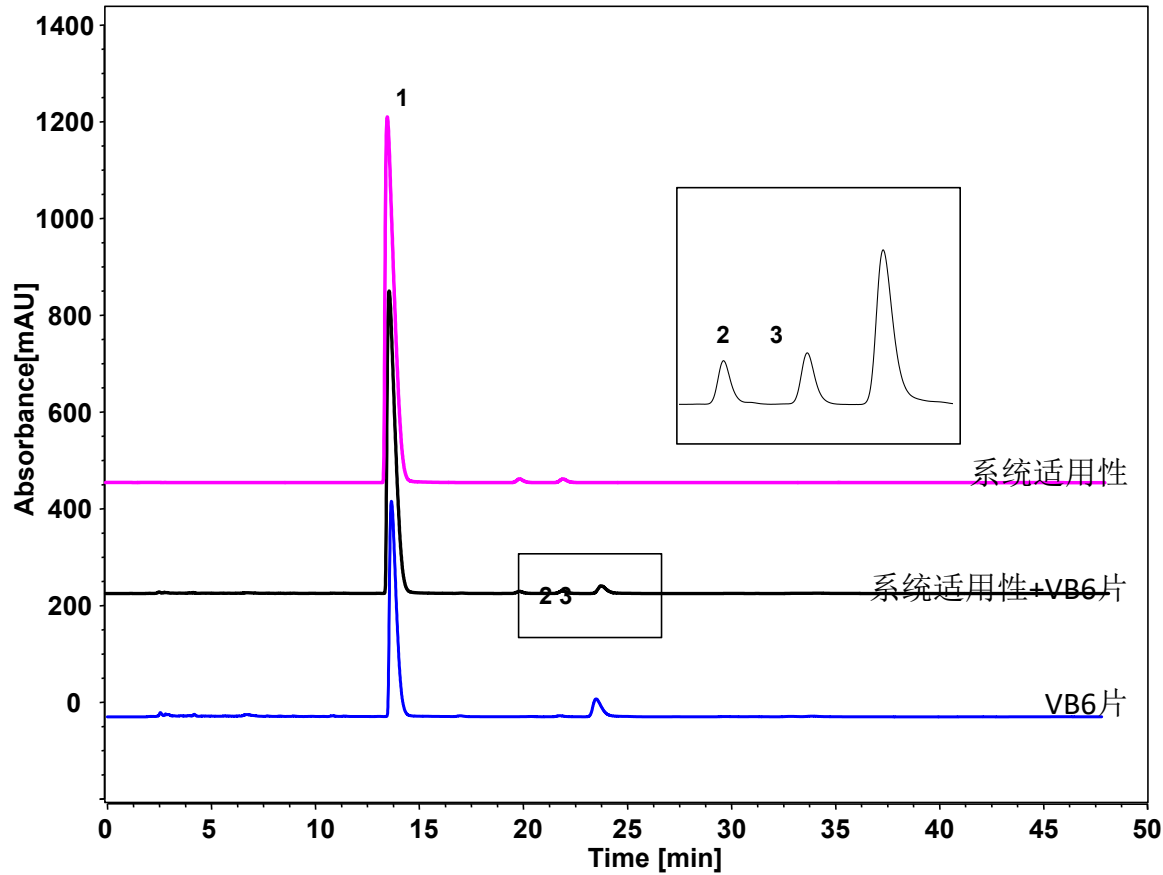