


# 左氧氟沙星 (中国药典)

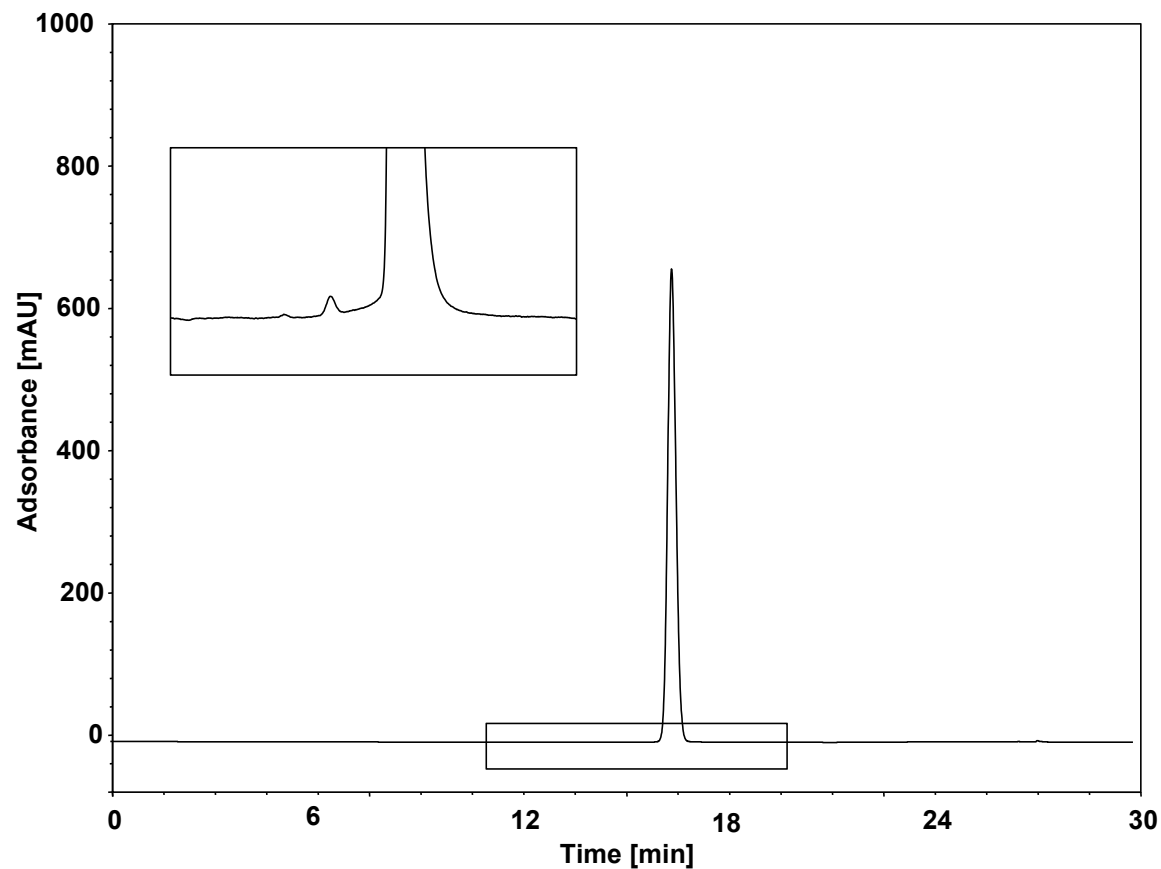

Column: ChromCore 120 C18, 3  $\mu$ m  
Dimension: 4.6  $\times$  250 mm  
Mobile Phase: A) 85/15 v/v 醋酸铵-高氯酸溶液, pH2.2/乙腈  
B) 乙腈  
Gradient: 

| t (min) | %A  | %B |
|---------|-----|----|
| -10     | 100 | 0  |
| 0       | 100 | 0  |
| 18      | 100 | 0  |
| 25      | 70  | 30 |
| 39      | 70  | 30 |
| 40      | 100 | 0  |

  
Flow Rate: 1.0 mL/min  
Temperature: 40  $^{\circ}$ C  
Injection: 10  $\mu$ L  
Detection: UV 254 nm  
Sample: 左氧氟沙星水溶液