




复方苯硝那敏片

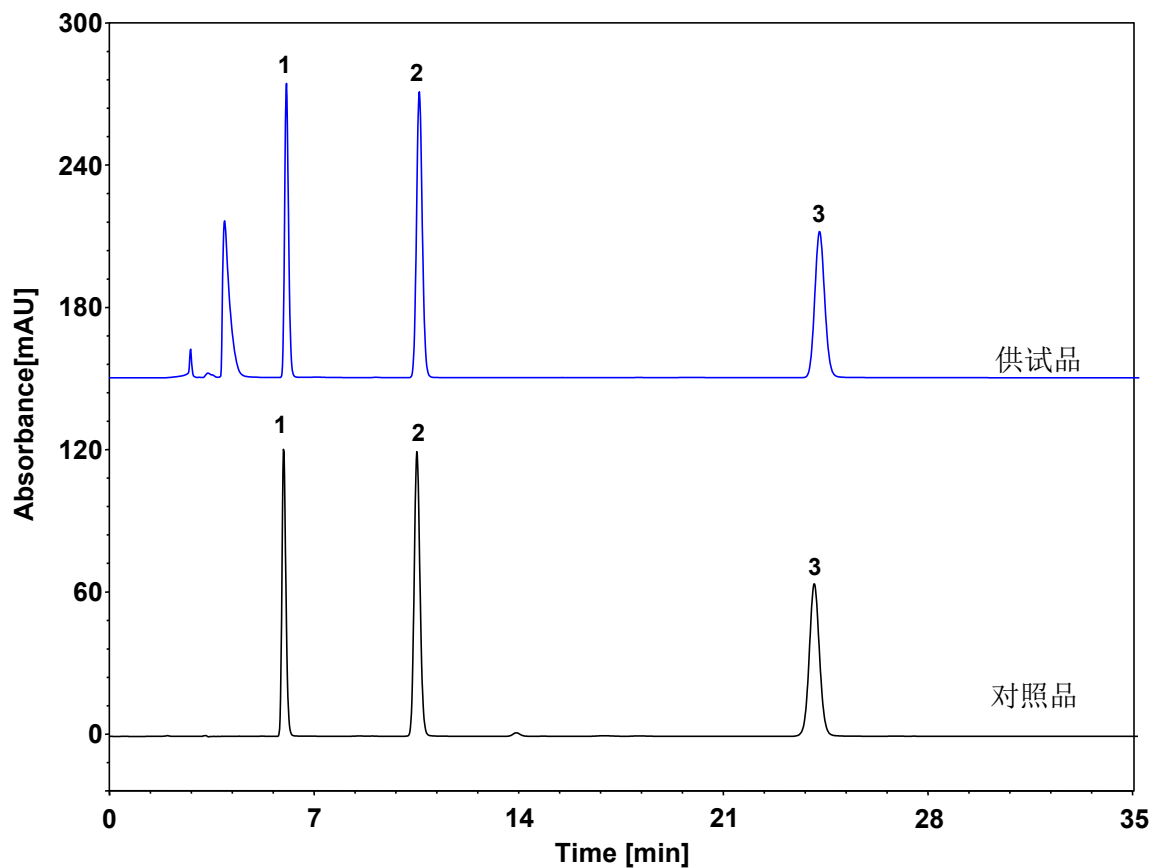


Column: ChromCore 120 C18, 5 μm
Dimension: 4.6 \times 250 mm
Mobile Phase: 60/40 v/v 甲醇/水
Flow Rate: 1.0 mL/min
Temperature: 25 $^{\circ}\text{C}$
Injection: 10 μL
Detection: UV 254 nm
Sample: 复方苯硝那敏片
Peaks:
1. 苯巴比妥
2. 硝西洋
3. 甲苯 (内标)

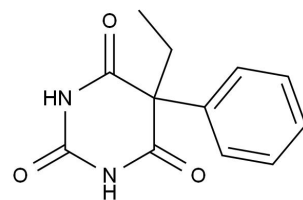

苯巴比妥

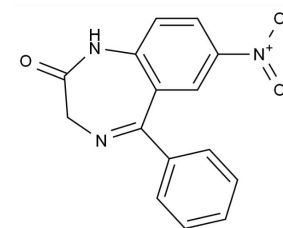

硝西洋

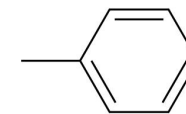

甲苯