




## 溴甲烷Bromomethane和卤化物Halogens在NanoChrom BP-624上的 分离与检测

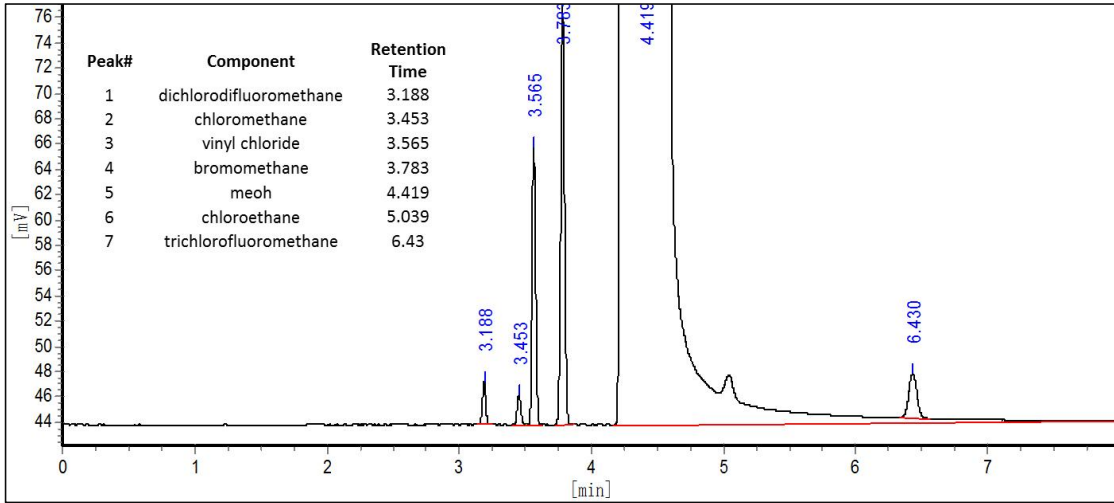

Column: NanoChrom BP-624  
 Dimension: 30m × 0.32mm × 1.80μm

Peak	Compound	Retention time (min)
1	dichlorodifluoromethane (二氯二氟甲烷)	3.188
2	chloromethane (氯甲烷)	3.453
3	vinylchloride (氯乙烯)	3.565
4	bromomethane (溴甲烷)	3.783
5	methanol (甲醇)	4.419
6	chloroethane (氯乙烷)	5.039
7	trichlorofluoromethane (三氯氟甲烷)	6.43