




二乙醇单异丙醇胺和三乙醇胺 在NanoChrom BP-INOWAX上的分离与检测

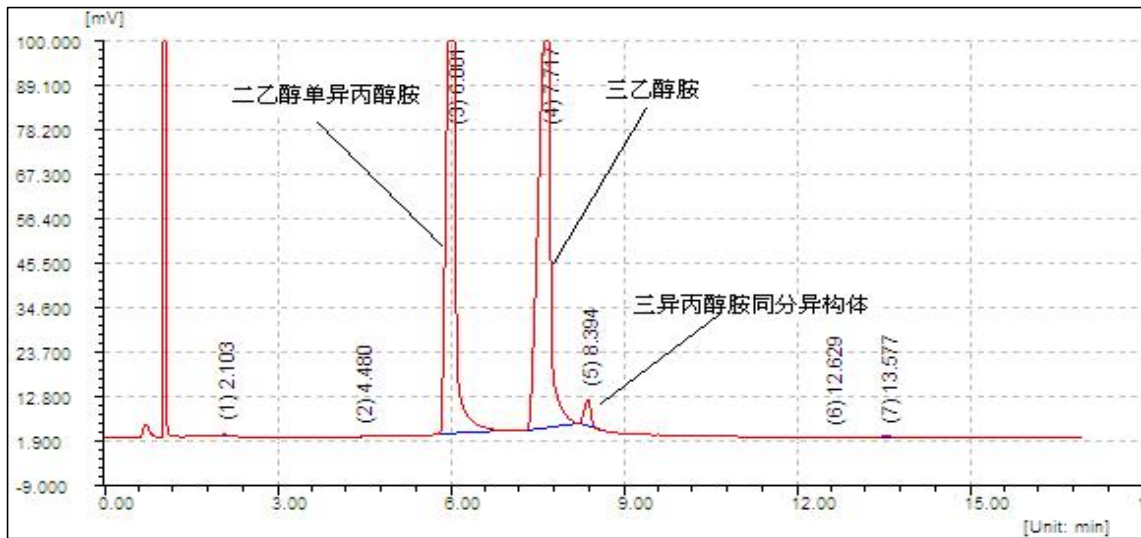

Column: NanoChrom BP-INOWAX
 Dimension: 30m×0.32mm×0.5μm
 Carrier: Nitrogen, head pressure: 0.14 MPa
 Inlet: Split, 280 °C, split flow: 30 mL/min
 Oven: 240 °C
 Inject: 1 μL
 Detector: FID 250 °C
 Peaks: 1. 二乙醇单异丙醇胺
 2. 三乙醇胺
 3. 三异丙醇胺同分异构体