




三氟氯菊酸/功夫酸在NanoChrom BP-FFAP上的分离与检测

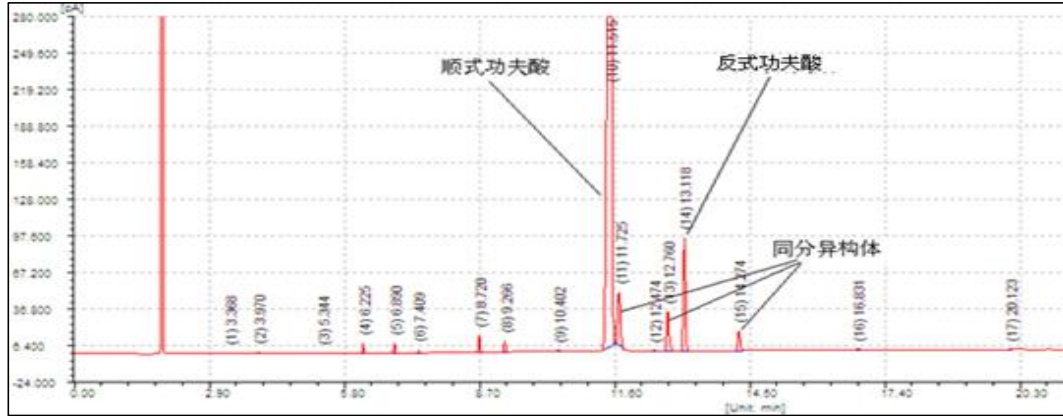

Column: NanoChrom BP-FFAP
 Dimension: 30m×0.32mm×0.25μm
 Carrier: Nitrogen, column flow: 1.5 mL/min
 Inlet: Split, 280 °C, split ratio: 80: 1
 Oven: 100 °C (hold 1 min) to 180 °C (hold 5.6 min) at 15 °C/min to 230 °C (hold 10 min) at 20 °C/min
 Inject: 0.4 μL
 Detector: FID 250 °C