


霉菌提取物

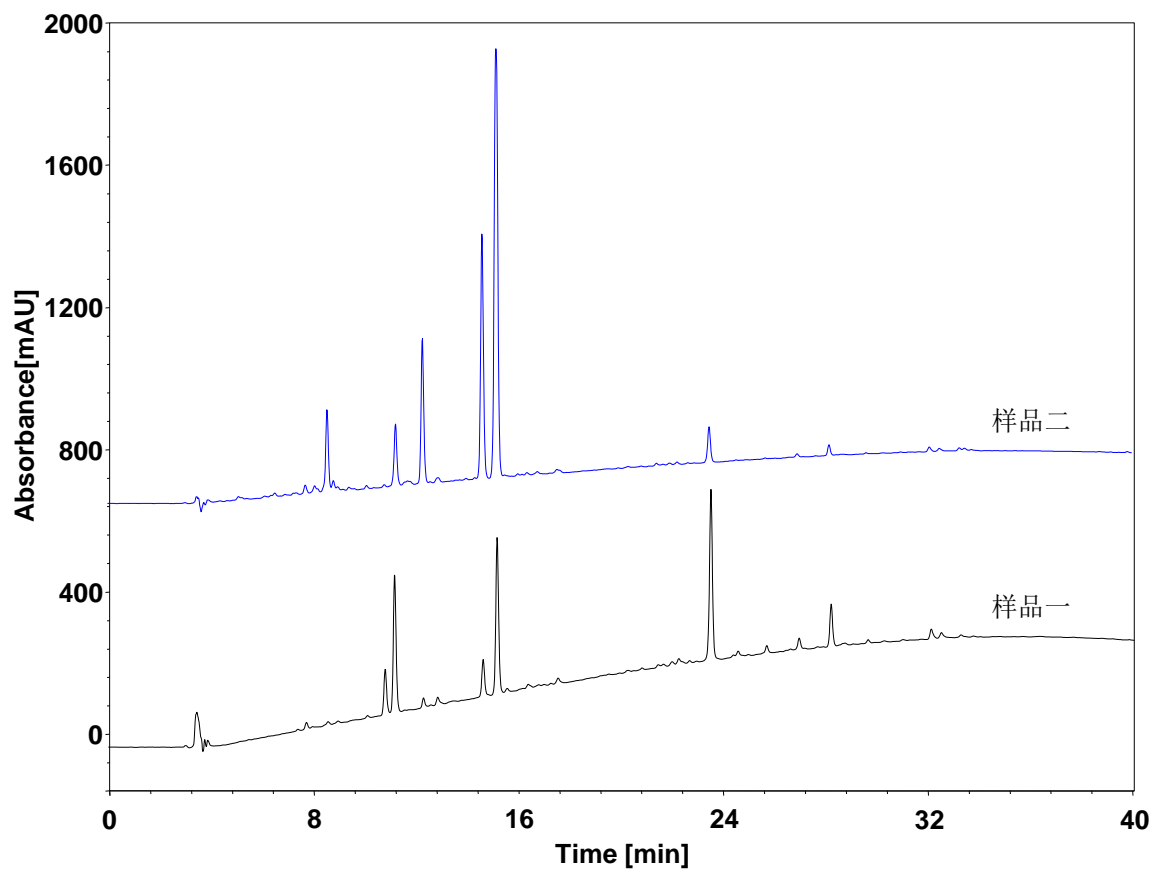

Column: ChromCore 300 C4-T, 5 μ m
Dimension: 4.6 \times 250 mm
Mobile Phase: A) 0.1% 三氟乙酸水
B) 0.1% 三氟乙酸乙腈

Gradient:	t (min)	%A	%B
	0	90	10
	40	5	95
	40.1	90	10
	55	90	10

Flow Rate: 1.0 mL/min
Temperature: 30 $^{\circ}$ C
Injection: 10 μ L
Detection: UV 214 nm
Instrument: Waters 2695
Sample: 霉菌提取物