


对乙酰氨基酚 (中国药典)

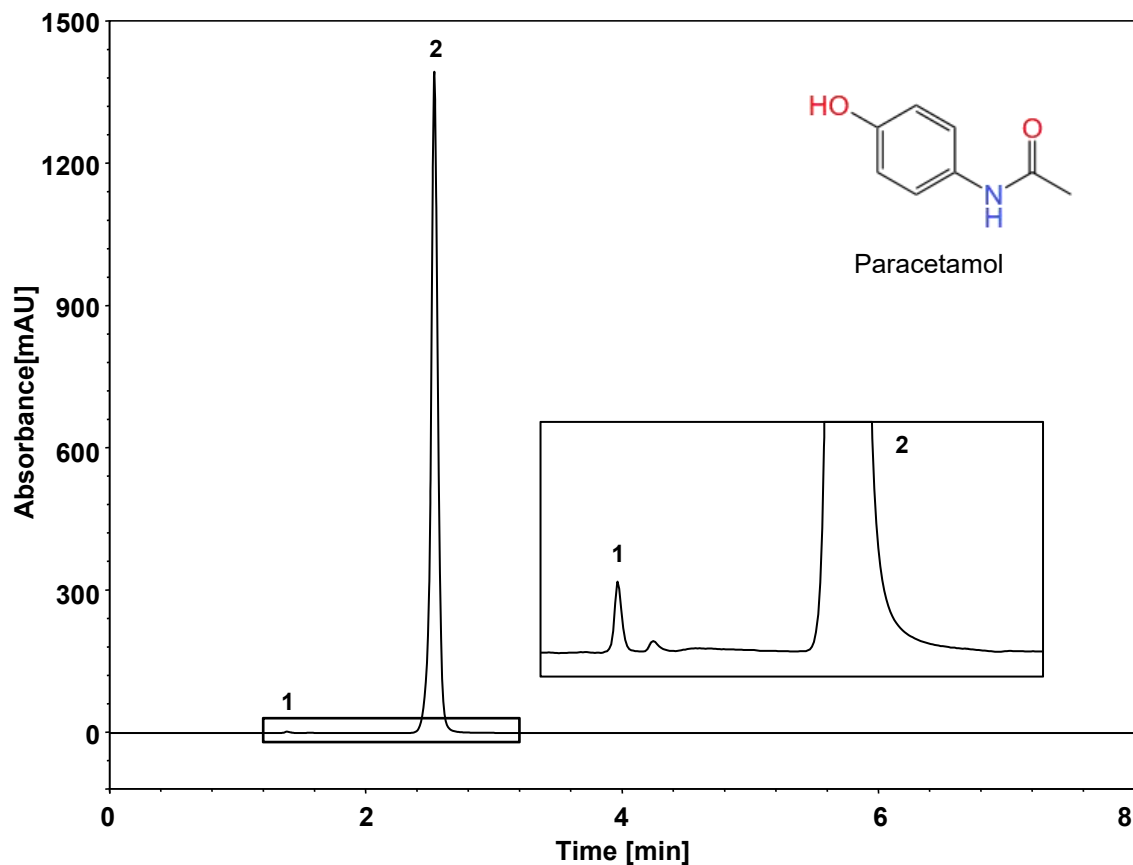

Column: ChromCore AQ C8, 1.8 μ m
Dimension: 2.1 \times 100 mm
Mobile Phase: 10/90 v/v MeOH/91mM phosphate buffer
Flow Rate: 0.3 mL/min
Temperature: 40 $^{\circ}$ C
Injection: 2 μ L
Detection: UV 245 nm
Instrument: Agilent Infinity II 1290 quaternary system
Peaks: 1. p-Aminophenol (对氨基苯酚)
2. Paracetamol (对乙酰氨基酚)

组分	保留时间 (min)	理论塔板数 (USP)	拖尾因子 (USP)	分离度 (USP)
p-Aminophenol	1.383	6926	1.13	/
Paracetamol	2.533	12204	0.95	14.54