


葛根 (中国药典)

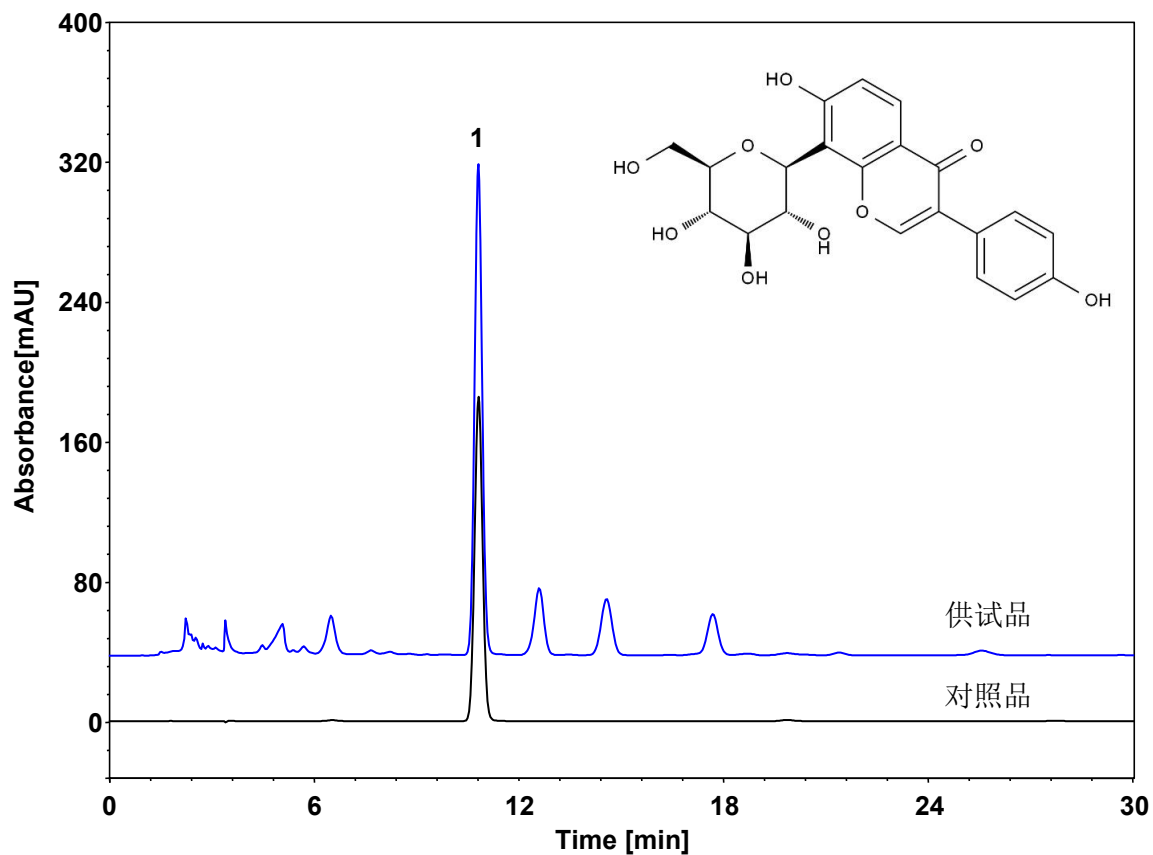

Column: ChromCore 120 C18, 5 μm
Dimension: 4.6 \times 250 mm
Mobile Phase: 25/75 v/v 甲醇/水
Flow Rate: 1 mL/min
Temperature: 40 $^{\circ}\text{C}$
Injection: 10 μL
Detection: UV 250 nm
Instrument: 岛津 LC-2030C
Sample: 葛根
Peak: 1. 葛根素

组分	保留时间 (min)	理论塔板数 (USP)	拖尾因子 (USP)	分离度 (USP)
葛根素	10.837	10808	1.04	/