


激素 (HCG)

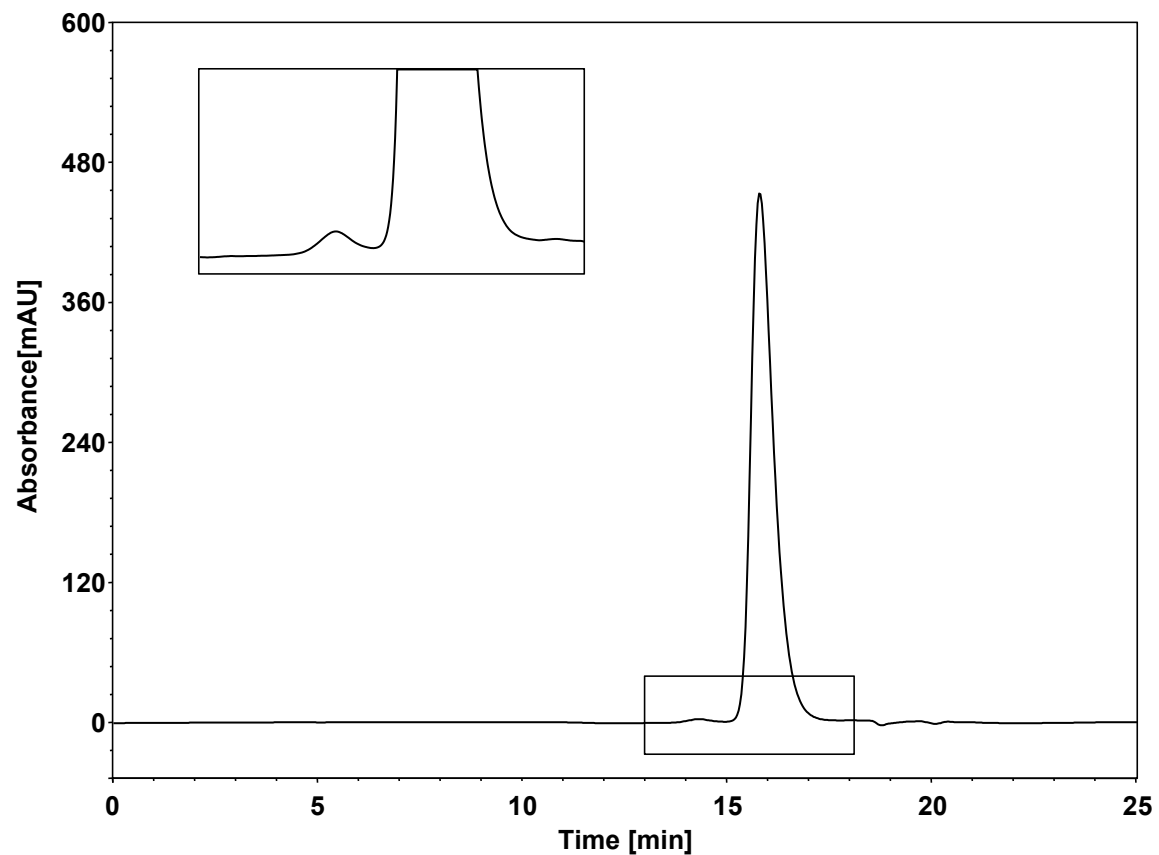

Column: BioCore SEC-300, 5 μ m
Dimension: 7.8 \times 300 mm
Mobile Phase: 10/90 v/v MeCN/ 200 mM Na₂SO₄ in 200 mM phosphate buffer, pH6.8
Flow Rate: 0.6 mL/min
Temperature: 25 $^{\circ}$ C
Injection: 25 μ L
Detection: UV 215 nm
Instrument: Waters 2695
Sample: 激素 (0.4 mg/mL MW 40 KD)