


盐酸达泊西汀

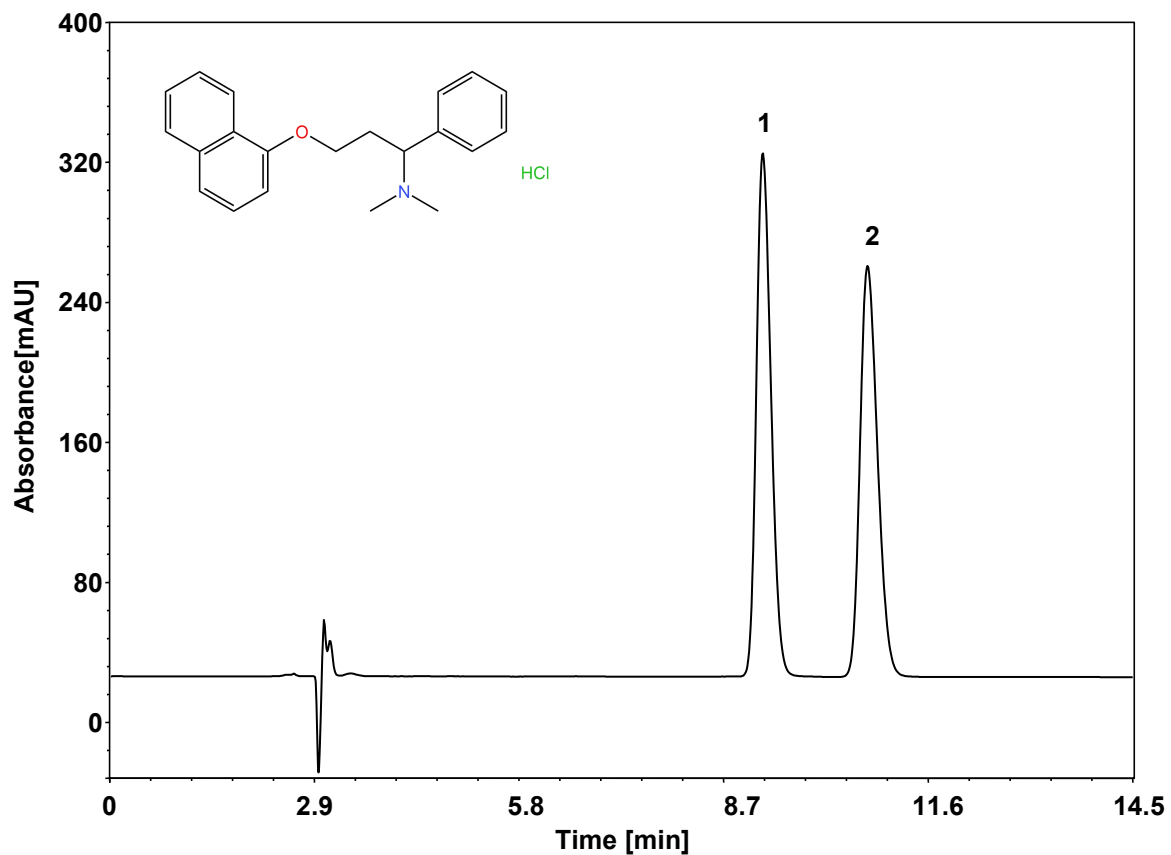

Column: ChiralCore Cel-J, 5 μ m
Dimension: 4.6 \times 250 mm
Mobile Phase: 90/10/0.1 v/v/v 甲醇/乙醇/二乙胺
Flow Rate: 1 mL/min
Temperature: RT
Injection: 10 μ L
Detection: UV 290 nm
Sample: 达泊西汀 (0.5 mg/mL, 乙醇溶解)

组分	保留时间 (min)	理论塔板数 (USP)	拖尾因子 (USP)	分离度 (USP)
1	9.249	9223	1.31	/
2	10.740	8626	1.30	3.52