


# 质粒

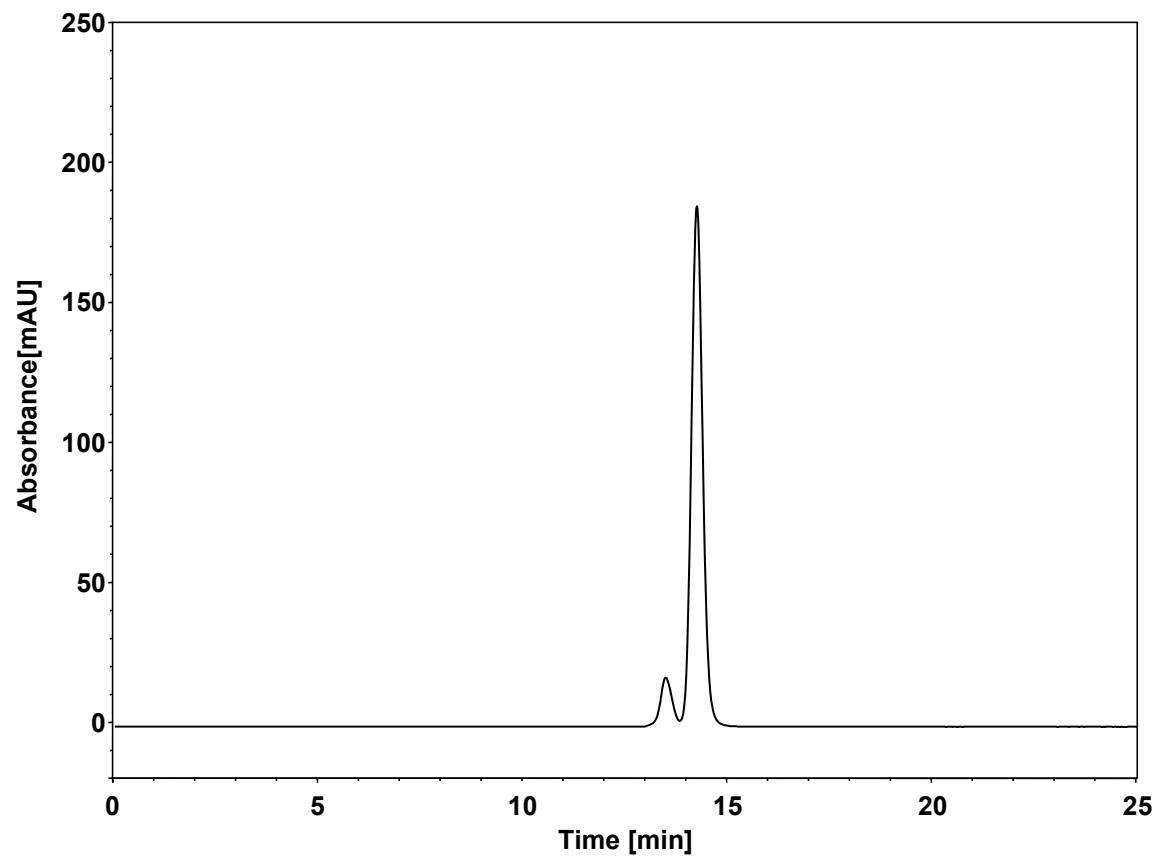

Column: BioCore SEC-1000, 5  $\mu\text{m}$   
Dimension: 7.8  $\times$  300 mm  
Mobile Phase: 300 mM NaCl in 50 mM phosphate buffer, pH6.8  
Flow Rate: 0.4 mL/min  
Temperature: 25  $^{\circ}\text{C}$   
Injection: 5  $\mu\text{L}$   
Detection: UV 260 nm  
Instrument: Agilent Infinity II 1260 (Bio)  
Sample: 质粒