


Labetalol

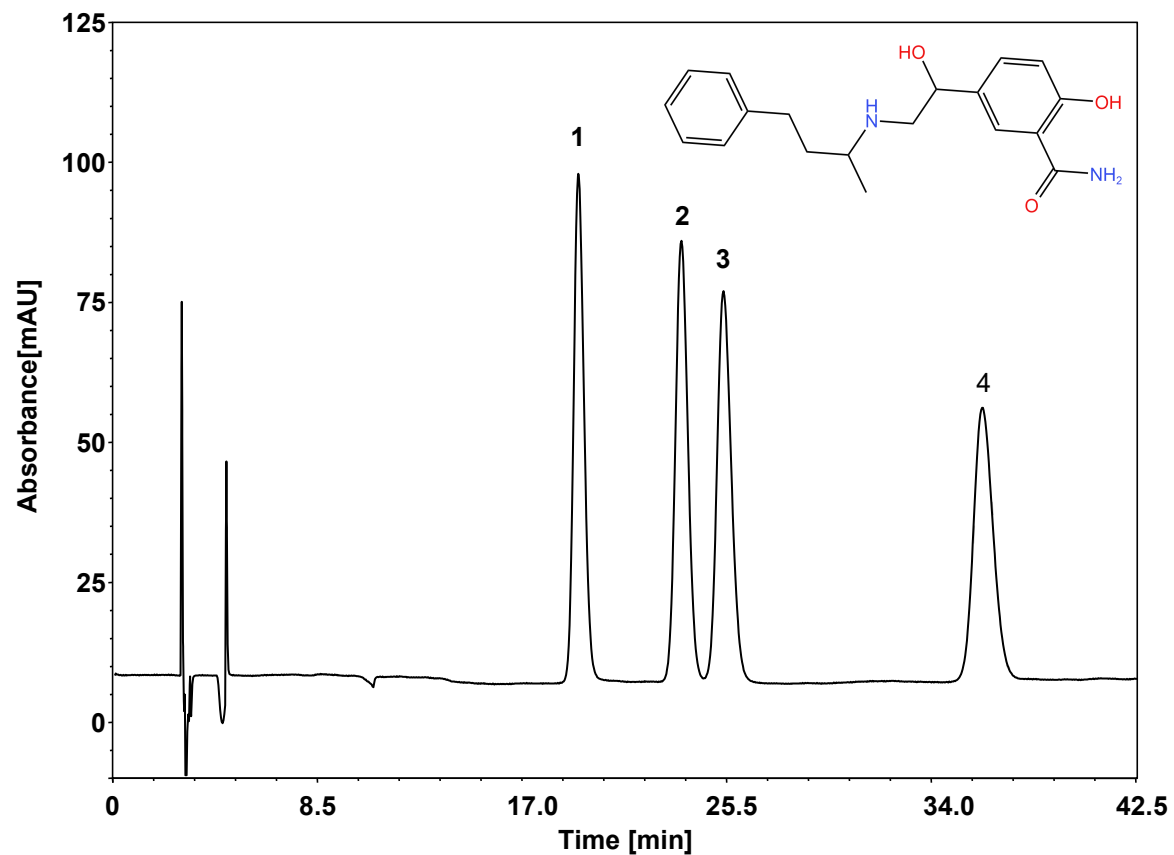

Column: ChiralCore Amy-D, 5 μm
Dimension: 4.6 \times 250 mm
Mobile Phase: 93/3/4/0.1 v/v 正己烷/乙醇/甲醇/乙醇胺
Flow Rate: 1.0 mL/min
Temperature: RT
Injection: 10 μL
Detection: UV 310 nm
Sample: Labetalol (2.5 mg/mL, 乙醇溶解)

组分	保留时间 (min)	理论塔板数 (USP)	拖尾因子 (USP)	分离度 (USP)
1	19.277	9270	1.29	/
2	23.560	9575	1.16	4.86
3	25.302	8652	1.13	1.70
4	36.077	9289	1.22	8.34